

TOPCOM®

EN | Instruction manual
NL | Gebruikershandleiding
FR | Manuel d'utilisation
DE | Bedienungsanleitung
ES | Manual de usuario
IT | Manuale utente
PT | Manual de utilizador
SV | Användarhandbok
TR | Kullanım kılavuzu
FI | Käyttöopas
NO | Brukermanual
PL | Instrukcja obsługi
CS | Návod na použití

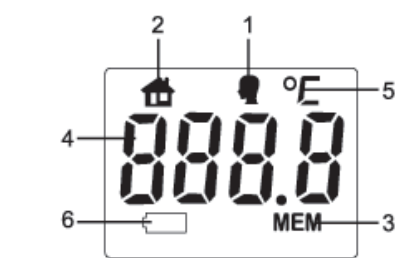TH-4676 | NON CONTACT INFRARED THERMOMETER

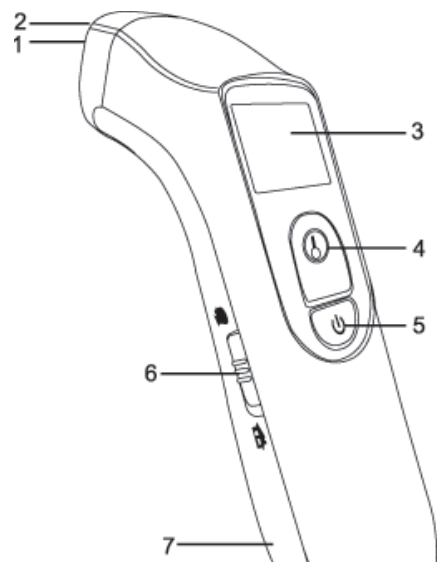smartwares
Europe

Jules Verneweg 87
5015 BH Tilburg
The Netherlands
www.smartwares.eu

www.smartwares.eu

EC REP

Shanghai International Trading Corp. GmbH (Hamburg)
Eiffestraße 80, 20537 Hamburg, Germany
Tel: 0049-40-2513175 Fax: 0049-40-255726
Email: shholding@hotmail.com

Shenzhen Aeon Technology Co., Ltd
RM6H02, Block 27-29
Tianxia Industrial Zone, Majialong
No. 133 of Yiyuan road, Nantou Street, Nanshan District
518052 Shenzhen | PEOPLE'S REPUBLIC OF CHINA

Storage & transport condition
Temperature: -25°C to 55°C. (-13°F to 131°F)
Storage humidity: 15% to 95% RH

EC REP	Authorized representative in the E.U.
	Manufacturer

EN Instruction manual

SAFETY INSTRUCTIONS

The following basic safety precautions should always be taken:

- Close supervision is necessary when the thermometer is used by, on, or near children, handicapped persons or invalids.
- Use the thermometer only for the intended use described in this manual.
- Do not use the thermometer if it is not working properly, or if it has suffered any damage.

PARTS DESCRIPTION

1. Infrared Sensor
2. Led Light
3. Display Screen
4. Measurement / Memory Button
5. ON/OFF Button
6. Mode Switch
7. Battery Cover

Display

1. Body mode indication
2. Object mode indication
3. Memory indication
4. Temperature reading
5. Temperature unit
6. Battery indication

BEFORE THE FIRST USE

- As with any thermometer, proper technique is crucial to obtaining accurate temperature readings. Please read this manual thoroughly and carefully before using.
- Always operate the thermometer in an operating temperature range from 10°C to 40° C and a relative humidity from 15% to 85%.
- Always store the thermometer in a cool and dry place (-20°C to 55°C) and a relative humidity from 15% to 90%.
- The device requires no calibration.
- The device contains no user serviceable parts.
- The user must check that the equipment functions safely and see that it is in proper working condition before use.
- No modification of this equipment is allowed.
- The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- Do not clean the device while it is being used.
- Avoid direct sunlight.
- Avoid dropping the thermometer, if it happens and you think the thermometer may be damaged, please contact customer service immediately.
- Do not touch the lens.
- Do not disassemble the thermometer.
- Basic safety precautions should always be observed, especially when the thermometer is used on or near children and disabled persons.
- This thermometer is not intended to substitute for a consultation with your physician
- This thermometer and the subject must remain in a stable environment for at least 30 minutes before measuring the temperature.
- Earwax in the ear canal or oil/grease on the forehead may cause an inaccurate measurement , so make sure the ear canal or forehead is clean.
- If you clean the ear or forehead, please wait for 5-10 minutes before measuring.
- When the measured temperature falls within the fever temperature range of ≥37.8°C (100.04°F) and <42.8°C (109.22°F), as indicated by the display, please consult with your physician immediately.
- Open the battery cover, remove the insulating plastic sheet and re-install the battery cover.
- The patient may be the operator of the thermometer. When the device is used to measure the temperature of a child, it has to be operated by an adult. Adults can measure their own temperature.
- The measured temperature is adjusted to show the body's core temperature.

Restrictions of use

- This thermometer is clinically proven to produce accurate temperature measurements. However, please be advised that the accuracy cannot be ensured when the thermometer is not clean. Check that the probe is clean before taking a measurement

INTENDED USE

- The infrared thermometer is intended for the intermittent measurement and monitoring of human body temperature from forehead measurement at home, clinics and hospital.
- A control measurement using a conventional thermometer is recommended in the following cases:
 - If the reading is surprisingly low.
 - For new-born infants up to 100 days old.
 - For children under three years of age who have a weakened immune system or who react unusually in the presence or absence of fever.

How does it work

The thermometer measures the infrared heat generated by the surface of the skin over the vessel and its surrounding tissue. The thermometer then converts it into a temperature value.

NOTE: The thermometer does not emit any infrared energy.

Display mode

- You can switch between two different modes with the mode switch:
- **Body mode:** This mode is used to measure the forehead temperature.
- **Object mode:** This mode is used to measure the object temperature.

Display color

Display color changes according to the measured temperature:

- **Green**(normal): < 37.5°C
- **Yellow**(mild fever): 37.5 – 37.9°C
- **Red** (high fever): > 38°C

Select the Temperature Unit

This meter provides two measurement units used for indicating the body/object temperature, °C or °F, for your preferred selection.

- Make sure the thermometer is OFF before selecting the temperature unit.

- Long-press the Measurement Button for 5 seconds until the °C/°F or °F/°C is displayed on the LCD panel as figure shown. Re-press the Measurement Button to select °C or °F unit.

Turn the Led Light ON or OFF

This meter provides Led Light to help users placing the thermometer at the correct position.

- Short-press the ON/OFF Button, the text "LED" or "OFF" will be displayed on the LCD panel.
- Quick press the measurement button, and select "LED or OFF" to turn the Led Light on or off.

Replacing the Battery

- The thermometer works with two 1.5 V AAA alkaline batteries. The thermometer will display an empty battery indication to alert you when the thermometer power is getting low, please follow the steps below to replace new batteries immediately.
- Remove the battery cover according to the arrow direction.
- Remove the old batteries and replace with two 1.5V AAA size alkaline batteries. Take care to match the Positive (+) and Negative (-) indications
- Close the battery cover according to the arrow direction.

NOTE:

- Although the thermometer works when the empty battery indication appears, we still recommend that you change the batteries to obtain an accurate result.
- Remove the batteries if stored for a long period of time.
- The batteries should be kept out of children's reach. If they are swallowed, promptly see a doctor for help.

USING THE DEVICE AS A BODY MEASUREMENT THERMOMETER

- Press the "ON/OFF" Button to turn on the thermometer first.
- Push the Mode switch to select Body mode. The device sounds one short beep when it is turned on and the temperature unit flashes to confirms the readiness to perform the measurement.
- Move the probe close to the forehead and take measurements. Make sure the probe is flat and close to the forehead, not at an angle. Perform a forehead measurement with a distance within 3 cm.
- Read the result. The measurement result will be done in 1 second. The reading is shown together with LED lighting, and one long beep informs about the temperature measurement, and after a while another, confirming the saving of the result to the memory and readiness for the next measurement.
- Press the "ON/OFF" button to turn off the unit, or leave it idle for 30 seconds, the unit will switch off automatically.

NOTE:

- As the forehead measurement temperature is likely to be affected by sweat, oil and the surroundings, the reading shall be taken as a reference only.
- If the probe is placed at an angle close to the forehead measurement, the reading will be affected by surrounding temperature. Babies' skin reacts very quickly to the ambient temperature. Therefore, do not take their temperature with the thermometer during or after breastfeeding, because the skin temperature may then be lower than the internal body temperature.
- If the measured temperature is < 32°C (89.6°F), the display will show with LO symbol.
- If the reading is <37.5°C (99.5°F) a green light will show.
- If the reading is between 37.5°C – 37.9°C (99.5°F – 100.2°F) a yellow light will show.
- If the reading is >37.9° (100.2°F) a red light will show.
- The thermometer will automatically turn off if left idle for 30 seconds.

As an Object Measurement Thermometer

- Press the "ON/OFF" Switch to turn on the thermometer first.
- Push the Mode switch to select Object mode. The machine emits two short beeps indicating that the object mode is turned on and the temperature unit flashes.
- Move the probe close to the object and take measurements. Make sure the probe is flat and close to the object, not at an angle. Perform a measurement within a distance of 3 cm. When ready, press the measurement button.
- Read the result. The measurement result will be done in 1 second. The reading is shown together with LED lighting , one long beep informs about the temperature measurement, and after a while another, short, confirming the saving of the result to the memory and readiness for the next measurement.
- Press the "ON/OFF" button to turn off the unit, or leave it idle for 30 seconds, the unit will switch off automatically.

Memory Function

- This thermometer stores the 25 most recent readings.
- Be sure the thermometer is OFF before recalling this memory.
- Press the Memory Button to enter the memory mode.
- Each time you press the Memory Button, a result will be displayed in the order of dates (latest result shown first), together with "MEM" and a number (from 1 to 25).
- The display will light up in green, yellow or red conforming to the memory reading.
- When the memory is full, the oldest result is deleted as the new one added. When the last record displayed in the display, press the Memory Button again to return to the first record.
- Press the ON/OFF Button to exit the memory function or leave the thermometer idle for 10 seconds to switch off automatically

Clearing the memory

- When in memory mode, long-press the ON/OFF Button for 3 seconds, release the button after the signal -- dEL displayed on the LCD panel as figure shown.
- Then long-press the ON/OFF Button for 3 seconds, release the button after the signal -- Cl displayed on the LCD panel as figure shown. All the memories would be deleted after that.
- **Silent mode:**
 - To measure temperature without 'beep' sound, please press and hold the measurement button for at least 7 seconds, until On or Off will display on the screen. By pressing the measurement button again you can switch between beep: on or off. Press the on/off button to confirm. When selecting 'off' the thermometer is in silent mode.

CLEANING AND MAINTENANCE

- The probe is not waterproof. Please wipe with a clean and dry cotton swab to clean the probe on the inside.
- The body of the thermometer is not water-resistant. Never put the thermometer under a running tap or submerge it into water. Use a soft and dry cloth to clean it. Do not use abrasive cleaners.
- Store the thermometer in a cool and dry location. Free from dust and away from direct sunlight.

ENVIRONMENT

- This appliance should not be put into the domestic garbage at the end of its durability, but must be offered at a central point for the recycling of electric and electronic domestic appliances. This symbol on the appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling of used domestic appliances you contribute an important push to the protection of our environment. Ask your local authorities for information regarding the point of recollection.

Support

You can find all available information and spare parts at www.smartwares.eu

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