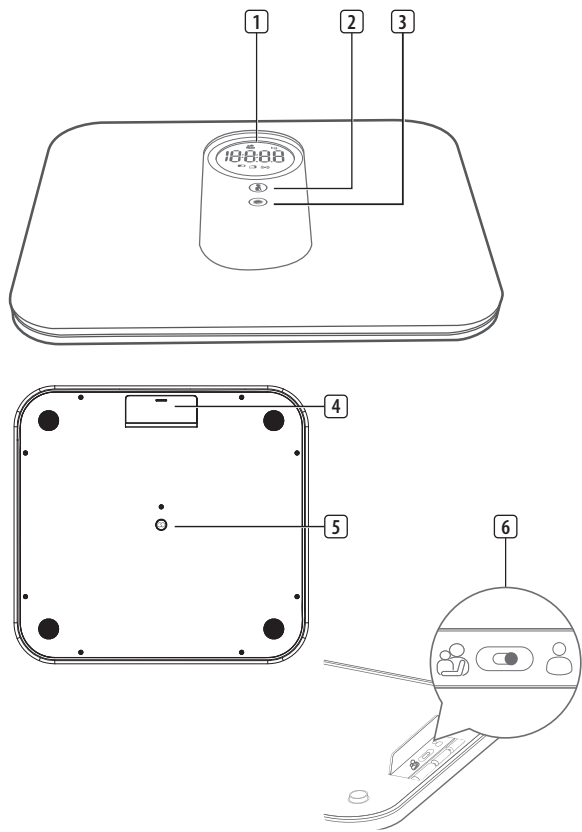
WG-2495

PARTS DESCRIPTION / ONDERDELENBESCHRIJVING / DESCRIPTION DES PIÈCES / TEILEBESCHREIBUNG / DESCRIPCION DE LAS PIEZAS / DESCRIZIONE DELLE PARTI / DESCRIÇÃO DOS COMPONENTES / BESKRIVNING AV DELAR / OPIS CZĘŚCI / POPIS SOUČÁSTÍ / POPIS SÚČASTÍ

[WEEE LOGO]

## EN Instruction manual

### SAFETY

- By ignoring the safety instructions the manufacturer can not be hold responsible for the damage.
- The appliance must be placed on a stable, level surface.
- Batteries**
  - Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
  - Keep batteries away from children. They are not a toy!
  - Do not open the batteries by force.
  - Avoid contact with metallic objects (Rings, nails, screws et cetera)there is a danger of short-circuiting!
  - As a result of a short-circuit batteries may heat up consider-ably or even catch fire. this may result in burns.
  - For your safety the battery poles should be covered with adhesive strips during transport.
  - If a battery leaks, do not rub your eyes with your hands. If in contact with the liquid, the hands should be washed thoroughly and eyes rinsed out with clear water, if symptoms persist please consult a doctor.

### PARTS DESCRIPTION

- LCD display
- Mother indicator
- Infant indicator
- Battery compartment
- Unit/ activation button
- Mother and infant mode switch

### BEFORE THE FIRST USE

- Take the appliance and accessories out the box. Remove the stickers, protective foil or plastic from the device.
- Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
- It is advisable to take measurement at the same time of the day and at the same place.
- Always remove your shoes and socks, it is advisable to discard as many as possible clothing, wipe clean your feet before using and dry them well. Never use the appliance with wet feet.
- Step on the weighing platform gently and when you are measuring, please keep still.
- The maximum weight capacity is 150kg. Graduation is 0,05kg.
- Remove the battery compartment cover at the back of device and place 3 x 1.5V AAA-size batteries. Replace the cover.
- Press the "UNIT" button at the back of the scale to select the measurement unit. The default measurement unit is "kg". You may press "UNIT" button to choose among kilogram, stone and pound.
- Install the App and Pair-Up Our App**
  - Note: It is only available for the device which supports Bluetooth 4.0 technology.
  - The App "Lifeseense international" could be freely downloaded at App Store and Google Play, you may search and install the app in your smart phone.
  - Turn on Bluetooth and App. Make sure both are ON when pair-up is proceeding.
  - Press and hold "UNIT" button in the back of the scale to start pair-up.

	Shown during pairing
	Shown if pairing succeed
	Shown if pairing failed

NOTE: a more detailed manual about how to use the "Lifeseense international" app can be found on: service.tristar.eu

#### Select Weighing Mode

The scale is equipped with two weighing modes: Normal Mode, Mother and Infant Mode. The weighing mode switch is located above the battery compartment. You may slide it to the left or the right to switch between Mother and Infant Mode and Normal Mode.

### USE

#### Normal Mode

- Slide the weighing mode switch to the right to enter Normal Mode.
- Place the scale flatly on a hard surface.
- Step onto the platform with barefoot and stand still until the measurement unit (e.g.kg)on the LCD stops blinking. The icon disappears and the icon blinks, the data will be transmitted automatically.
- The scale will display the measuring result for about ten seconds and then shut off automatically. If it fails, the icon lights, the measurement data will be saved in the scale.

	this icon will blink during transmitting, the data will be transmitted automatically, if the data is complete the icon will disappear
	If succeed the symbol disappears, if failed the symbol remains

#### Mother and infant mode

- Slide the weighing mode switch to the left to enter Mother and Infant Mode.
- Place the scale flatly on a hard surface.
- Step onto the platform with barefoot. The scale powers on automatically and the symbol Mother and Infant Mode lights up.
- Then the mother indicator lights up. Please stand still until the measurement unit (e.g. kg) on the LCD stops blinking. The measurement data will be transmitted to your smartphone automatically.
- About 3 seconds after your weight is locked, both the mother indicator and the infant indicator light up. The scale restores to zero automatically.
- Hold the baby in your arms and step onto the platform. Please stand still until the measurement unit (e.g. kg) on the LCD stops blinking.
- When the baby's weight is locked, the mother indicator goes out while the infant indicator keeps on. Your baby's weight will be transmitted to the smartphone and displayed on the LCD for about 10 seconds.
- Then the scale powers off automatically.

#### Error prompt

Error	Description	Solution
	Overload, the device will power off in four seconds.	Stop using this scale for measurement.
	Low Battery. The device will power off in four seconds.	Replace all three batteries at the same time. Please purchase the authorized batteries for replacement.
	Failure of pairing up. or failure of activation.	Please check: Bluetooth is ON, App is ON, both devices are within the transmission distance of Bluetooth.

### CLEANING AND MAINTENANCE

- Do not immerse the appliance in water or any other liquid.
- Clean the appliance with a damp cloth. Never use harsh and abrasive cleaners, scouring pad or steel wool, which damages the device.

### GUARANTEE

- This product is guaranteed for 24 months granted. Your warranty is valid if the product is used in accordance to the instructions and for the purpose for which it was created. In addition, the original purchase (invoice, sales slip or receipt) is to be submitted with the date of purchase, the name of the retailer and the item number of the product.
- For the detailed warranty conditions, please refer to our service website: www.service.tristar.eu

### ENVIRONMENT

Fixed width narrow. This appliance should not be put into the domestic garbage at the end of its durability, but must be offered at a central point for the recycling of electric and electronic domestic appliances. This symbol on the appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling of used domestic appliances you contribute an important push to the protection of our environment. Ask your local authorities for information regarding the point of recollection.

### Support

You can find all available information and spare parts at service.tristar.eu!

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