
- EN | Instruction manual
- NL | Gebruikershandleiding
- FR | Manuel d'utilisation
- DE | Bedienungsanleitung
- ES | Manual de usuario
- IT | Manuale utente
- PT | Manual de utilizador
- SV | Användarhandbok
- TR | Kullanım kılavuzu
- FI | Käyttöopas
- NO | Brukermanual
- PL | Instrukcja obsługi
- CS | Návod na použití

## TH-4656 | EAR & FOREHEAD THERMOMETER

AUGUST 2019
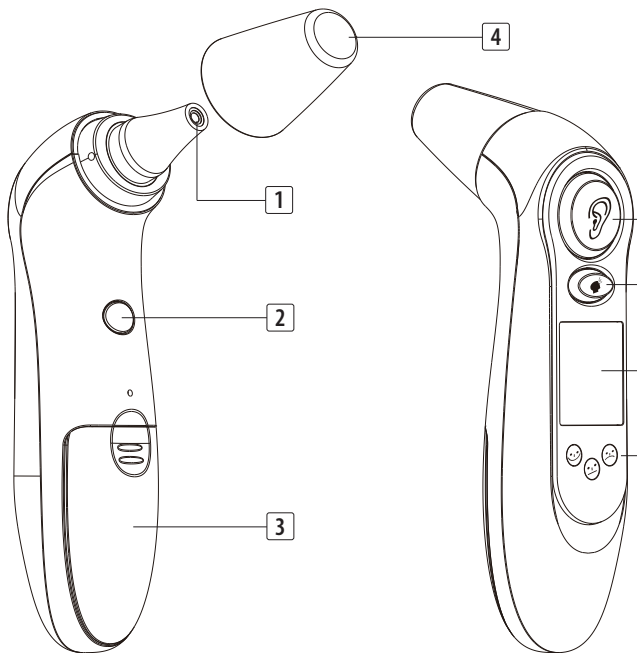
## EN Instruction manual

### SAFETY INSTRUCTIONS

The following basic safety precautions should always be taken:

- Close supervision is necessary when the thermometer is used by, on, or near children, handicapped persons or invalids.
- Use the thermometer only for the intended use described in this manual.
- Do not use the thermometer if it is not working properly, or if it has suffered any damage.

### PARTS DESCRIPTION

- |                          |                              |
|--------------------------|------------------------------|
| 1. Infrared sensor       | 5. Ear measuring button      |
| 2. Memory and Set button | 6. Forehead measuring button |
| 3. Battery cover         | 7. LCD display               |
| 4. Protection cap        | 8. LED indicator             |

### BEFORE THE FIRST USE

- As with any thermometer, proper technique is crucial to obtaining accurate temperature readings. Please read this manual thoroughly and carefully before using.
- Always operate the thermometer in an operating temperature range from 10°C to 40°C and a relative humidity from 15% to 95%.
- Always store the thermometer in a cool and dry place (-25°C to 55°C) and a relative humidity from 15% to 95%.
- The device requires no calibration.
- The device contains no user serviceable parts.
- The user must check that the equipment functions safely and see that it is in proper working condition before use.
- No modification of this equipment is allowed.
- The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- Do not clean the device while it is being used.
- Avoid direct sunlight.
- Avoid dropping the thermometer, if it happens and you think the thermometer may be damaged, please contact customer service immediately.
- Do not touch the lens.
- Do not disassemble the thermometer.
- Basic safety precautions should always be observed, especially when the thermometer is used on or near children and disabled persons.
- This thermometer is not intended to substitute for a consultation with your physician.
- This thermometer and the subject must remain in a stable environment for at least 30 minutes before measuring the temperature.
- When the measured temperature falls within the fever temperature range of ≥37.8°C (100.4°F) and <42.9°C (109.2°F), as indicated by the red LED on display, please consult with your physician immediately.
- Dry the device in a ventilated area if water enters the case from the outside. After moisture has evaporated, you can resume use.
- When the operating environment is relatively dry, strong electromagnetic interference usually occurs. When this happens, the device may be affected as follows:
  - The device turns off;
  - The device restarts;
  - The above phenomenon does not affect the basic safety and essential performance of the device, and the user can use it according to the instruction. If you want to avoid the above phenomenon, please use it according to the environment specified in the manual.

### Restrictions of use

- This thermometer is clinically proven to produce accurate temperature measurements. However, please be advised that the accuracy cannot be ensured when the thermometer is not clean. Check that the probe is clean before taking a measurement

### INTENDED USE

The infrared thermometer is intended for the intermittent measurement and monitoring of human body temperature at home.. A control measurement using a conventional thermometer is recommended in the following cases:

- If the reading is surprisingly low.
- For new-born infants up to 100 days old.
- For children under three years of age who have a weakened immune system or who react unusually in the presence or absence of fever.

### Temperature indication

There are three temperature indication-LED's in this unit.The relevant LED will light up when the measured human body temperature is between 34.0°C(93.2°F) and 42.9°C(109.3°F).

- Green LED: 34.0°C(93.2°F)~ 37.4°C(99.4°F)
- Orange LED: 37.5°C(99.5°F)~37.9°C(100.3°F)
- Red LED: 38.0°C(100.4°F)~42.9°C(109.3°F)

### How to measure forehead temperature

1. Remove protective cap before use. Place the thermometer in your hand with your thumb on the Forehead button.
2. Hold the thermometer so that the probe is aimed at the center of the forehead, midway between eyebrow and hairline. The probe doesn't touch the forehead, keep the distance within1cm.
3. Press the Forehead button the start the measurement. Hold the device steady until you hear two beeps. You can now read the temperature on the display.
4. You can take a new measurement after 2 seconds.
5. Press the MEM/SET button to return to standby mode. The thermometer will automatically go into standby mode after 1 minute.

### How to measure ear temperature

1. Remove protective cap before use. Place the thermometer in your hand with your thumb on the Ear button.
2. Insert the probe very gently and slowly into the ear canal. Point the probe toward the eardrum inside the ear canal.
3. Press the Ear button to start the measurement. Hold the device steady until you hear two beeps. You can now read the temperature on the display.
4. You can take a new measurement after 6 seconds. Always take measurements in the same ear, to avoid discrepancies.

### How to measure object temperature

1. Remove protective cap before use. Place the thermometer in your hand with your thumb on the Ear and Forehead buttons.
2. Short-Press the Ear and Forehead button simultaneously until "food" is shown on the display.
3. Hold the thermometer so that the probe is at a distance of 1 cm of the object whose temperature you want to measure.
4. Press the Ear button to start the measurement. Hold the device steady until you hear two beeps. You can now read the temperature on the display.
5. Press the Ear and Forehead button simultaneously again to return to standby mode..

### Checking the memory

There are a total of 9 memory storages for recording temperature measurements. The current measurement is always stored in the last storage place. When all storages are occupied, the oldest measurement will be overwritten.

- Press the MEM/SET button to check the memory data. Press it multiple times to check older measurements.
- Press the Ear and Forehead buttons at the same time to return to standby mode.

### Setting the date and time

- Press and hold the MEM/SET button for 2 seconds to open the settings menu.
- Use the Forehead button to switch between the following settings:
  - 12hrs (AM/PM) or 24hrs
  - Hour
  - Minute
  - Year
  - Month
  - Day
  - Sleep mode on/off
- Adjust the setting by pressing the MEM/SET button multiple times, and confirm with the Forehead button.

### Switching between Celsius and Fahrenheit

- When in standby or in measuring mode, press and hold the Ear and Forehead buttons until the temperature switches from Celsius to Fahrenheit, or vice versa.

### Replacing the Battery

- When the low battery indicator appears on the LCD , or the thermometer does not respond at all, you should replace the battery as soon as possible.
- Slide off the battery cover, remove the old battery, and replace it with a new CR2032 battery. Make sure the polarity off the battery is correct when placing it.
- **NOTE:**
  - Although the thermometer works when the empty battery indication appears, we still recommend that you change the batteries to obtain an accurate result.
  - Remove the batteries if stored for a long period of time.
  - The batteries should be kept out of children's reach. If they are swallowed, promptly see a doctor for help.
  - When the thermometer displays one of the following error messages: "ErrP" or "ErrH", this could indicate a hardware problem. Please contact the manufacturer if this occurs.
  - When the thermometer displays the following error message: "ErrE", please check if the operating temperature is within the allowed parameters. Wait for 10 seconds and take another measurement.

### CLEANING AND MAINTENANCE

- To ensure accurate measurements it is very important to keep the probe tip clean and free of scratches.
- Finger prints, earwax or dirt will affect the accuracy of the thermometer. In order to get an accurate measurement, please clean the probe tip.
- Gently wipe its surface with a cotton swab slightly moistened with alcohol and immediately wipe dry with a clean cotton swab.
- Use a soft dry cloth to clean the thermometer. Never clean the thermometer with an abrasive cleaner, thinner, benzene, or submerge the thermometer into water or other liquids.

### ENVIRONMENT

This appliance should not be put into the domestic garbage at the end of its durability, but must be offered at a central point for the recycling of electric and electronic domestic appliances. This symbol on the appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling of used domestic appliances you contribute an important push to the protection of our environment. Ask your local authorities for information regarding the point of recollection.

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
TH-4656 | EAR & FOREHEAD THERMOMETER  
AUGUST 2019

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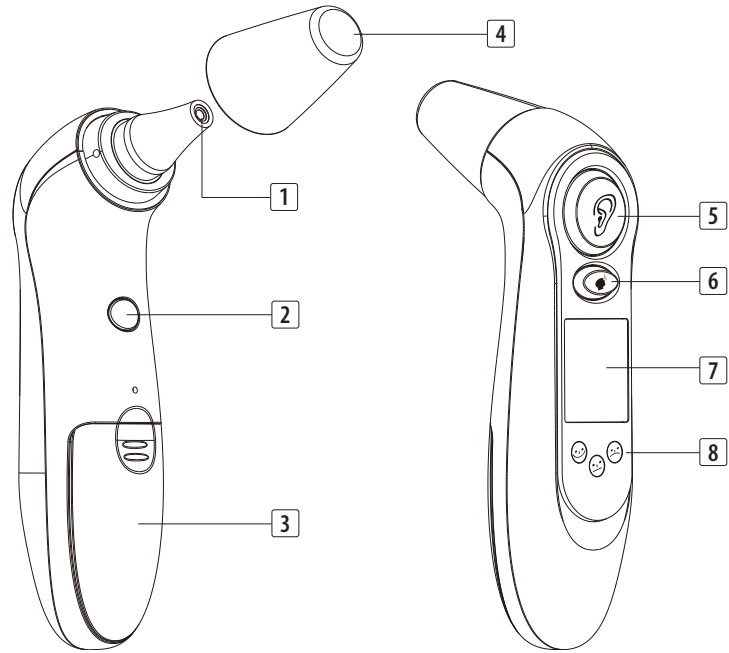**EMC compliance**  
Declaration - electromagnetic emission  
The device is intended for use in the electromagnetic environment specified below. The customer or the user of device should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The device is suitable for domestic establishment and in establishment directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	

Declaration - electromagnetic immunity  
The device is intended for use in the electromagnetic environment specified below. The customer or the user of device should assure that it is used in such an environment.

Immunity test	IEC 60601test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tie. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/ burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0.5kV, ± 1 kV line(s) to lines ± 0.5kV, ± 1 kV, ± 2 kV line(s) to earth	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270°and 315°  0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0°  0 % UT; 250/300 cycles	Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: UT is the ac mains voltage prior to application of the test level.  
Declaration - electromagnetic immunity  
The device is intended for use in the electromagnetic environment specified below. The customer or the user of device should assure that it is used in such an environment.

Immunity test	IEC 60601test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3V 0.15 MHz to 80MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 MHz	Not applicable	Portable and mobile RF communications equipment should be used no closer to any part of device, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance
Radiated RF IEC 61000-4-3	10V/m 80 MHz to 2.7 GHz	10V/m	$d = 1.2\sqrt{P}$ 150 KHz to 80 MHz $d = 1.2\sqrt{P}$ 80 MHz to 800 MHz $d = 2.3\sqrt{P}$ 80 MHz to 2.7 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range.  Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.  
NOTE 2 These guidelines may not apply in all situations.  
Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.  
Field strengths from fixed RF transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the

location in which device is used exceeds the applicable RF compliance level above, device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating device.  
Over the frequency range 0.15 MHz to 80 MHz, field strengths should be less than 3 V/m.  
Recommended separation distances between portable and mobile RF communications equipment and device  
The device is intended for use in an electromagnetic environment in which radiated RF disturbance are controlled. The customer or the user of device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and device, as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m	0.15 MHz to 80 MHz	80 MHz to 800 MHz	80 MHz to 2.7 GHz
		$d = 1.2\sqrt{P}$	$d = 1.2\sqrt{P}$	$d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23	
0.1	0.38	0.38	0.73	
1	1.2	1.2	2.3	
10	3.8	3.8	7.3	
100	12	12	23	

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer  
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.  
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

smartwares Europe	NL: 088 5940501 BE: 078 481094(NL) FR: 0235 560 650	Lokaal tarief Lokaal tarief Tarif local
Jules Verneweg 87 5015 BH Tilburg The Netherlands	DE: +49 (0) 1805 010762 14Ct./Min. dt. Festnetz, Mobil max. 40Ct./Min.	
www.topcom.eu	UK: +44 (0) 345 230 1231 ES: +34 938427589	

Disposal in accordance with Directive 2012/19/EU(WEEE).  
 Refer to instruction manual

**SN** Serial number.  
 Applied apart of type BF

**CE** 0197  
Complies with the European Medical Device Directive (93/42/EEC) and amended by directive 2007/47/EC requirements. Notified body TÜV Rheinland (0197).

**IP22**  
The first number 2: Protected against solid foreign objects of 12,5mm φ and greater. The second number: Protected against vertically falling water drops when enclosure tilted up to 15°. Vertically falling drops shall have no harmful effect when the enclosure is tilted at any angle up to 15°, on either side of the vertical.

Shenzhen Dongdixin Technology Co., LTD.  
Floor 1-2, No.3 Building, Fanshen Xusheng Industrial Estate, Xilixiaobaomang, Nanshan District, Shenzhen 518108, China  
Tel: 0086-755-27652316  
E-mail: service@e-caretalk.com  
www.e-caretalk.com

**EC REP**  
Shanghai International Holding Corp. GmbH (Europe)  
Eiffestraße 80, 20537 Hamburg, Germany  
Tel: 0049-40-2513175 Fax: 0049-40-255726

<b>EC REP</b>	Authorized representative in the E.U.
	Manufacturer

When the operating environment is relatively dry, strong electromagnetic interference usually occurs. At this time, the device may be affected as follows:  
- the device stops output;  
- the device turns off;  
- the device restarts;  
The above phenomenon does not affect the basic safety and essential performance of the device, and the user can use it according to the instruction. If you want to avoid the above phenomenon, please use it according to the environment specified in the manual.